


OMNI TH

Tissue Homogenizer

- Ergonomic handheld design or can be post mounted
- Uses Omni Tip Plastic probes or Stainless Steel Probes
- Variable motor speeds from 5,000 to 35,000 rpm
- Processing range from 0.03ml-2L

The **Omni Tissue Homogenizer (TH)** is ergonomically designed to reduce hand strain during extended periods of operation, or it can be post-mounted to an optional stand assembly. The Omni TH processes most samples in less than 30 seconds.

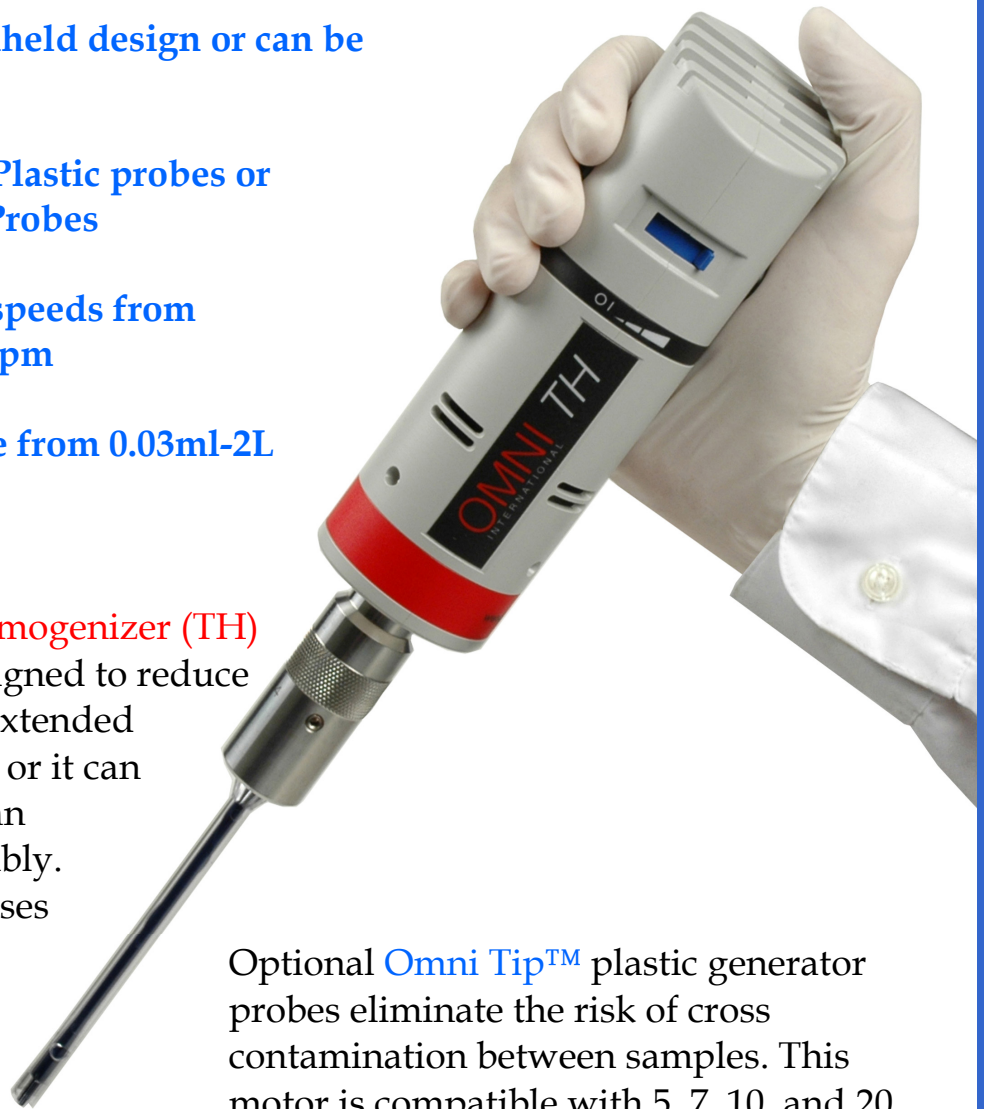

Optional **Omni Tip™** plastic generator probes eliminate the risk of cross contamination between samples. This motor is compatible with 5, 7, 10, and 20 diameter generator probes.



OMNI
International
The Homogenizer Company

1000 Williams Drive, Suite 1024
Marietta, GA 30066
omni@omni-inc.com

800-776-4431
770-421-0206 Fax
www.OMNI-INC.com

OMNI TH

Tissue Homogenizer

Features

- Powerful 125 watt motor
- Light weight ergonomic design
- Quiet operation, under 72db
- Variable speed control from 5,000 to 35,000 rpm
- 1/4 turn bayonet probe mounting system for fast release
- CSA and CE approved

Specifications

Processing Range: 0.03ml-2L

Speed: Variable from 500 to 35,000 rpm

Power Rating: 125 watt Max

Dimensions: 16cm (6.5 in) height, 5.1cm (2 in) diameter

Weight: 538g (19 oz)

Sound Level: <68 db

Warranty: 1 year motor warranty

Standards Approval: CSA and CE approved



Omni TH shown with optional motor stand assembly and stainless steel probe.

The **Omni Tissue Homogenizer** operates with **Omni Tip™** disposable/reusable plastic probes or stainless steel generator probes.



OMNI
International
The Homogenizer Company

1000 Williams Drive, Suite 1024
Marietta, GA 30066
omni@omni-inc.com

800-776-4431
770-421-0206 Fax
www.OMNI-INC.com